FARETEC INC.
RESPONDER PRODUCTS

EMERGENCY CARE | MILITARY | CIVILIAN

# RESPONDER PELVIC BINDER






**STABILIZES**
PELVIC FRACTURES

**OPTIMAL FORCE**
TENSION ADJUSTMENT

**RAPID APPLICATION**
SAVES TIME

1x **SINGLE PATIENT USE**
DISPOSABLE

SCAN FOR
**TRAINING VIDEO**
faretec.com/ifu

**TRUSTED BY PROFESSIONALS**

**DESIGNED FOR RELIABILITY**

**BUILT TO PERFORM**

FDA | CE/MDR | ISO13485
CAGE: 0R123
Made in U.S.A.

IFU-EN-V2.0
MAY 2024

# RESPONDER PELVIC BINDER

## INSTRUCTIONS FOR USE



### INTENDED USE

The Responder Pelvic Binder is a non-invasive, circumferential pelvic belt intended to stabilize pelvic fractures during transport to a definitive care facility.

### INDICATIONS FOR USE

The Responder Pelvic Binder is indicated for application by prehospital and medical professionals or trained personnel, to individuals who have, or are suspected of having, a pelvic ring fracture where hemorrhage might occur.

### CONTRAINDICATIONS

Do not utilize the Responder Pelvic Binder if anatomical landmarks are absent.

## WARNINGS

- Use only for the stabilization of suspected or confirmed pelvic ring fractures.
- Improper application or over-tightening may cause additional injury or patient discomfort.
- Do not delay transport to a definitive care facility to achieve perfect reduction.
- Monitor the patient's condition continuously. This device does not stop internal bleeding.
- Remove the binder if it interferes with life-saving procedures.

## PRECAUTIONS

- Remove objects from patient's pockets and clothing prior to application.
- Monitor the patient's skin for pressure injury while in use.
- 1x Device is single use only.
- Reuse of device may cause illness (e.g., infection) or injury (reduced performance).
- Dispose of used device per healthcare facility biohazard protocol.

## CLINICAL INFORMATION FOR USERS

Familiarization and training on the Responder Pelvic Binder, adherence to healthcare organization evidence-based standards of care and clinical guidelines including use of PPE, are required for safe use of this product. Failure to utilize this device in a manner consistent with these instructions for use, training, and within clinical best practice guidelines may result in serious illness or injury.


# RESPONDER PELVIC BINDER

## PRODUCT ANATOMY

### BUILT FOR PERFORMANCE

Engineered to rapidly stabilize pelvic fractures with adjustable, circumferential compression.

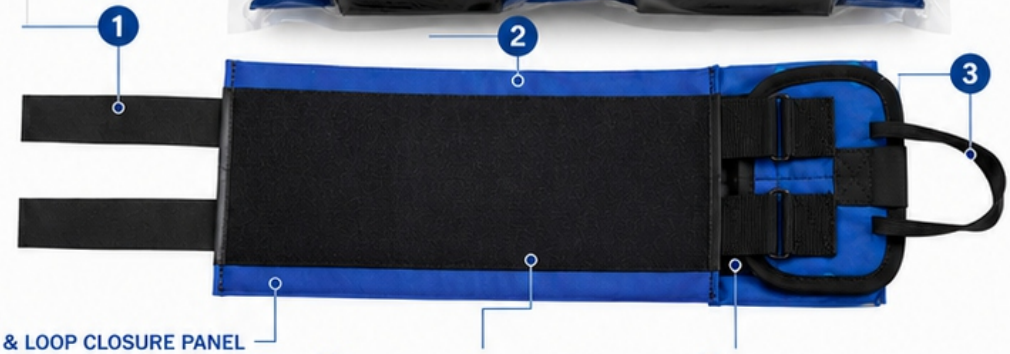


**1 COMPRESSION STRAPS**
Dual adjustable straps apply circumferential compression to stabilize the pelvis.

**2 MAIN BINDER**
Durable, wide wrap design distributes force evenly and comfortably.

**3 COUNTERFORCE HANDLE**
Provides counterforce during tensioning to achieve secure and controlled compression.

**4 HOOK & LOOP CLOSURE PANEL**
Large closure area ensures strong, secure attachment and easy adjustment.

**5 INTEGRATED PULL LOOPS**
Reinforced loops for quick grip and efficient strap tensioning.

**6 SIZING ADJUSTMENT AREA**
Excess material can be folded or trimmed to achieve the proper fit for each patient.

### INTENDED USE

The **Responder Pelvic Binder** is a non-invasive, circumferential pelvic belt intended to stabilize pelvic fractures during transport to a definitive care facility.

### INDICATIONS FOR USE

The Responder Pelvic Binder is indicated for application by prehospital and medical professionals or trained personnel, to individuals who have, or are suspected of having, a pelvic ring fracture where hemorrhage might occur.

### CONTRAINDICATIONS

Do not utilize the Responder Pelvic Binder if anatomical landmarks are absent.

### WARNINGS

- Use only for the stabilization of suspected or confirmed pelvic ring fractures.
- Improper application or over-tightening may cause additional injury or patient discomfort.
- Do not delay transport to a definitive care facility to achieve perfect reduction.
- Monitor the patient's condition continuously. This device does not stop internal bleeding.
- Remove the binder if it interferes with life-saving procedures.

### PRECAUTIONS


- Remove objects from patient's pockets and clothing prior to application.
- Monitor the patient's skin for pressure injury while in use.
- Device is single use only.
- Reuse of device may cause illness (e.g., infection) or injury (reduced performance).
- Dispose of used device per healthcare facility biohazard protocol.

FARETEC INC
RESPONDER PRODUCTS



## DIMENSIONS:

SIZE:
LENGTH - 53”in (134 cm)
WIDTH - 8”in (20 cm)

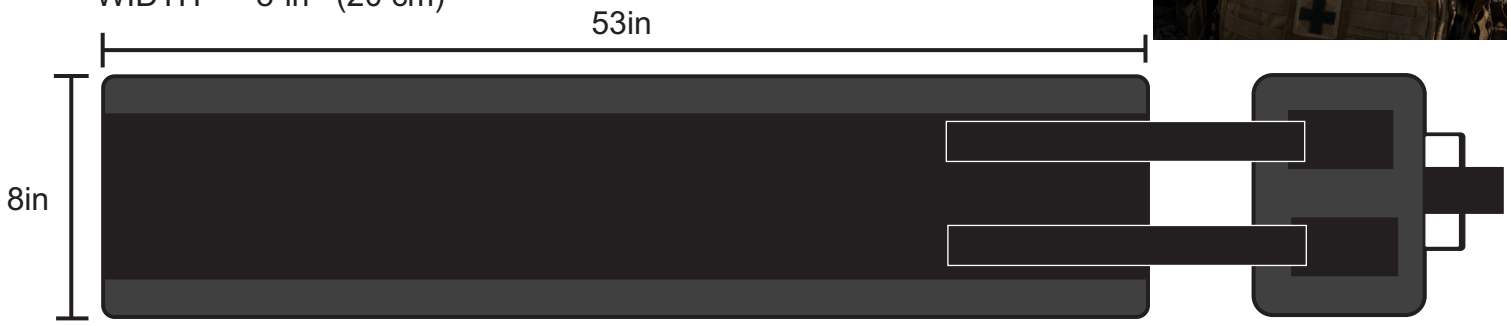



OUTSIDE

STRAP LENGTH: 20”in (50cm)



INSIDE

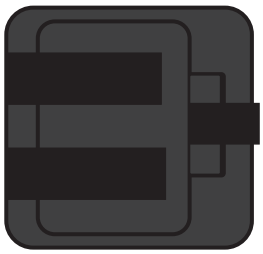

PACKED DIMS:
- 9 x 8 x 1.3” inches
- FOLDED & VACUUM SEALED

WEIGHT: 1 lb, .45 kg
MATERIAL: CORDURA
MilSpec Velcro, 2” NYLON

Once removed from packaging, the binder is ready for immediate use. Just unroll and slid under the patient as per the IFU.

